



INSTRUCTION SHEET

HOW TO TAKE MEASUREMENTS OF A VELOCE WHEELCHAIR

Material(s) required:

- None

Tool(s) required:

- Measuring tape
- Digital protractor

1. **SEATH WIDHT** : From the outside of the frame minus $\frac{1}{4}$ ".




2. **SEAT DEPTH** : From the front of the back cane to the front of the seat upholstery.



3. **FRAME LENGTH** : From the front of the back cane to the beginning of the frame guard.



4. **SEAT TO BACK ANGLE** : From the seat frame to the back cane.

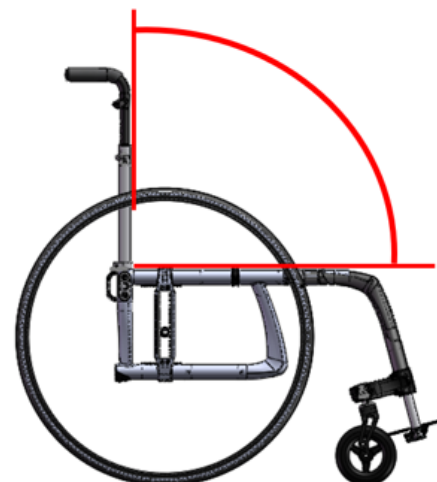

5. FOOTREST ANGLE :



6. FOOTPLATE ANGLE :



7. SEAT TO FOOTPLATE LENGTH : From the front of the seat upholstery to the footplate support tube.




8. ARMREST HEIGHT : From the top of the seat frame to the top of the armrest pad.



9. BACK CANE HEIGHT :


A. All straight and 8° backrest canes (aluminium & carbon) :

STANDARD : From the bottom of the clamp to the start of the bend of the push handle.

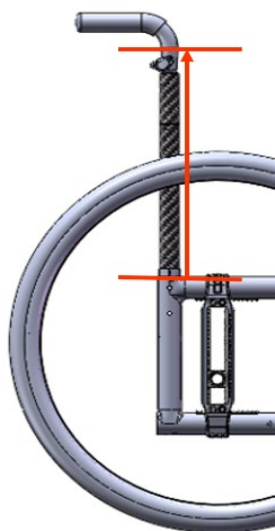

B. All backrest cane fixed straight and 8° (aluminium & carbon) :

FIXED : From the bottom of the clamp to the top of the cap.

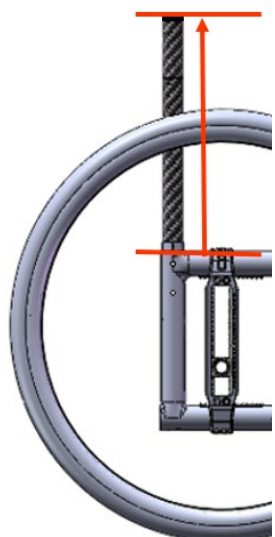

C. All angle adjustable backrest canes right angle and 8° (aluminium & carbon) : ANGLE ADJUSTABLE : From the top of the clamp to the bend of the push handle.

