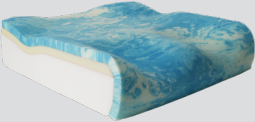

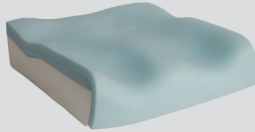



# SEAT CUSHION OPTIONS WITH CLINICAL JUSTIFICATION



PRODUCT	FOAM	WEIGHT CAPACITY	CLINICAL JUSTIFICATION
<b>BIOFIT</b> 	<b>3 LAYERS:</b> <ul style="list-style-type: none"> <li>Bottom: firm molded layer with SIT technology*</li> <li>Middle: medium density visco foam</li> <li>Top: soft density smartGEL infused visco foam</li> </ul> 	<ul style="list-style-type: none"> <li>• 300 lbs (136.1 kg)</li> <li>• 22-24" (55.8-61 cm) widths up to 400 lbs (181.4 kg)</li> </ul>	<ul style="list-style-type: none"> <li>• Client requires moderate skin protection &amp; positioning.</li> <li>• Top layer smartGEL infused visco technology reduces heat buildup, can minimize skin integrity issues.</li> <li>• Contours include notches for greater trochanters, sacral relief and pre-ischial bar to provide proper pelvic positioning to assist with trunk positioning for MRADL's and improved sitting tolerance throughout the day.</li> <li>• NXT SIT (Selective Immersion Technology) perforations support ischial/sacral area to prevent sliding and shear, as well as providing immersion for bony prominences.</li> <li>• Beveled front for foot propellers or to accommodate tight hamstrings.</li> </ul>
<b>NÜFIT</b> 	<b>3 LAYERS:</b> <ul style="list-style-type: none"> <li>Bottom: firm molded layer with SIT technology*</li> <li>Middle: medium density visco foam</li> <li>Top: soft density foam for comfort</li> </ul> 	<ul style="list-style-type: none"> <li>• 300 lbs (136.1 kg)</li> </ul>	<ul style="list-style-type: none"> <li>• Client requires mild-moderate skin protection &amp; positioning.</li> <li>• Contours include notches for greater trochanters, sacral relief and pre-ischial bar to provide proper pelvic positioning to assist with trunk positioning for MRADL's and improved sitting tolerance throughout the day.</li> <li>• NXT SIT (Selective Immersion Technology) perforations support ischial/sacral area to prevent sliding and shear, as well as providing immersion for bony prominences.</li> <li>• Beveled front for foot propellers or to accommodate tight hamstrings.</li> </ul>

\*See [website](#) for complete description of SIT technology. \*\*See [website](#) for more information on Smartx3D.

## NOTES

---



---



---



---



---

SCAN HERE FOR  
SIZE RANGES  
AND FURTHER  
PRODUCT DETAILS.



# COVER OPTIONS



FOR SEAT CUSHIONS	MATERIAL	CLINICAL JUSTIFICATION
COVER Smartx3D (STANDARD COVER)	• 4-way stretch, Smartx3D "spacer fabric"	<ul style="list-style-type: none"><li>• Enhanced immersion</li><li>• Breathable</li><li>• Moisture wicking</li><li>• Vibration damping</li><li>• Three layer technology creates breathable channel and "microclimate"</li><li>• The channel between layers provides cushioned feel for increased comfort</li></ul>
COVER INCONTINENCE	• Hydrotuff and startex	<ul style="list-style-type: none"><li>• Protects foam from fluids</li><li>• Wipeable for ease of cleaning</li></ul>

NOTES

---

---

---

---

---

---

---

---

---

---